NOOANCE

PARIS

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 **NOOANCE INNOVATIONS**
22 Rue Beaujon, 75008 Paris, France

Table of Contents

The Science Of Light Therapy	28
Company Information	29
Product Information	30
Device Specifications	
Device Accessories	
Adverse Reaction	
Contraindications	
Warnings	
Precautions	
Travel And International Use	
Disclaimer	
Risks	38
Potential Risks	
Light Sensitivity Test	39
How To Use The Device (Two ways to use)	
Connect the controller to use	
Connect the adapter to use	
Installation of Straps	
Treatment Protocol	
Charging	
Care And Storage	
Cleaning	
Troubleshooting	
Disposal	
Warranty	45
Symbols	46

The Science Of Light Therapy

The skin has the wonderful ability to absorb light energy, the skin utilizes this energy to kick-start the cells into producing more cellular fuel (ATP), boosting their performance and allowing them to do what they do naturally, but at an enhanced rate.

Super-charged cells from exposure to light therapy use this energy to prompt skin fibroblasts to produce increased collagen, elastin proteins and stimulate the body's natural healing and rejuvenating processes.

Light therapy is recognized medically and also proven clinically, and offers many therapeutic benefits that are safe, non-invasive and chemical free.

Company Information

Company

NOOANCE INNOVATIONS

Address

22 Rue Beaujon, 75008 Paris, France

Product Information


Device Specifications

LED wavelengths	633 nm +/- 2 nm (Red) 830 nm (Infrared)
Max. output optical power density	2 x 45 mW/cm2
Operating temperature	0°C to 40°C 32°F to 104°F
Storing and transport temperature	-5°C to 55°C 23°F to 131°F
Operating humidity	30% to 95% without condensation
Storage and Transport Humidity	20% to 95% without condensation
Charging supply	5V = 2A
Charging port	USB-C
Charging cable	USB-C to USB-A
Charging cable length	2 meters / 6.562 feet
Charging adapter power output	2A
Battery	3.7V, 3600mAh
Battery type	Lithium-ion
Charging time	240 minutes / 4 hours
Recommended use	1 time per day
Treatment time	10 minutes
Dimensions of the device	331*276.7*5 mm
Weight of the device	0.712 kg / 25.12 oz

The red light of the LEDs may appear orange after a few seconds of use as the eye adapts to this brightness. This is completely normal.

Device Accessories

1. NOOANCE Led Mask Multi-Functions
2. Manual
3. Adjustable Straps
4. Power Supply
5. EU Adapter
6. US Adapter
7. UK Adapter
8. USB Cable
9. Controller
10. Storage Bag

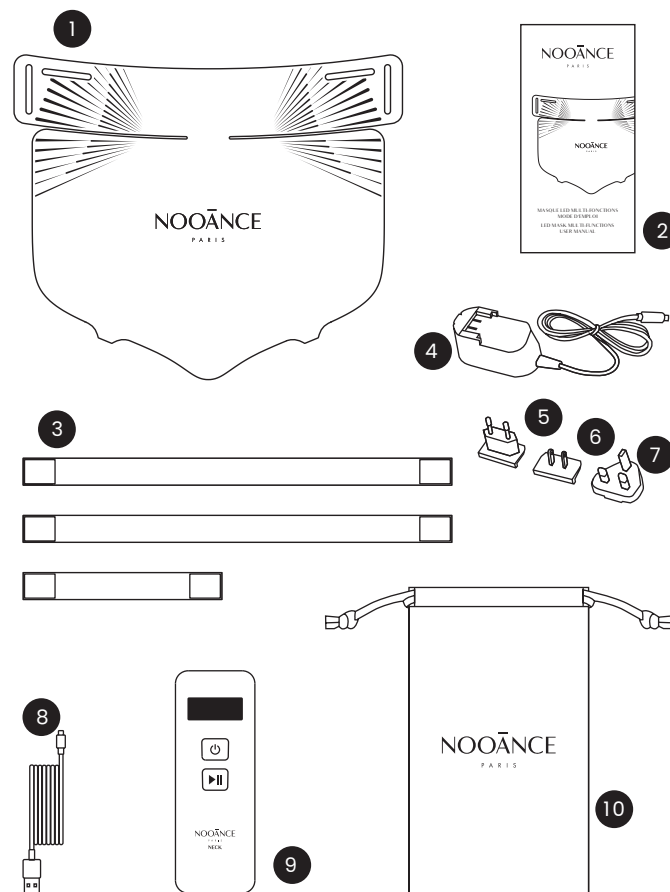Adverse Reaction

If you have any adverse reaction, such as discomfort, skin discoloration, skin irritation or other concerns, from using this device, stop using it immediately and consult your health professional.

Skin irritation, bruising, contact burns, dizziness, nausea, and headaches are rare but possible side effects. Immediately stop using the device and contact your General Practitioner.

The red light of the LEDs may appear orange after a few seconds of use as the eye adapts to this brightness. This is completely normal.

Contraindications

- DO NOT use the device if you have any active cancer or cancerous lesion.
- DO NOT use the device if you are pregnant, planning on becoming pregnant, or breastfeeding.
- DO NOT use the device on people under the age of 18.
- DO NOT use the device over areas with phlebitis, thrombophlebitis, and varicose veins.
- DO NOT use the device if you have a known allergic reaction to silicon. The device contains plastics and silicon that come into contact with the skin during use. When in doubt of allergies first contact your General Practitioner before using the device.
- DO NOT use the device if you suffer from any photosensitive disorder (sensitization to light).
- DO NOT use the device if you are taking any medication or drugs that can cause photosensitivity. If you are not sure about your medication contact your health professional.
- DO NOT use the device if you suffer from light-induced headaches.
- DO NOT use the device if you have a medical history of seizures or epilepsy that was triggered by light.

Note: For reasons of hygiene, we do not recommend sharing your device with anyone else.

Note: We highly suggest taking caution when taking any (skin related) damaging medications. If in doubt, consult your General Practitioner.

Note: If you are under the care of a physician, consult your physician before using this device.

⚠ Warnings

- STOP using the device and consult your General Practitioner if you experience an adverse reaction from using the device.
- STOP using the device if you experience symptoms of hypotension (dizziness or light-headedness) while using the device.
- STOP using the device if it is overheating. Turn off the device immediately.
- STOP your treatment if it feels hot or your skin feels uncomfortably warm.
- DO NOT let children play with the device. Cleaning and charging the device shall not be done by children.
- DO NOT use the device to treat rosacea, moles, warts, open sores, cancerous lesions, or any skin condition.
- DO NOT use the device on spots without pigmentation or areas on the skin that have vitiligo.
- DO NOT shine the light from the LEDs directly into your eyes or someone else's eyes. This may cause temporary loss of sight.
- DO NOT place the device if you have any open wounds.
- DO NOT use the device outdoors to prevent damage to the product.
- DO NOT sit, lay down or step on the device. This will break the device.
- DO NOT drop or insert any object into any openings of the device.
- DO NOT use the device if it is not working properly, or if it has been dropped or damaged. For warranty and repairs contact customer support.
- DO NOT use or store the device near hot surfaces.
- DO NOT charge the device when it is close to or might get in contact with water.
- DO NOT submerge the device in water.

- DO NOT place the device under direct sunlight.
- DO NOT attempt to modify the device in any way. Modifying or opening the device will void your warranty. The device contains no serviceable parts. This device contains batteries that are non-replaceable.
- DO NOT perform maintenance or take apart the product. In case of malfunctioning or damage, contact customer support.
- DO NOT use the device if damaged or when the charging cable is damaged.
- DO NOT use the device with non-manufacturer-approved accessories.
- DO NOT use solvents or abrasive cleaners on the device.
- DO NOT use the device in any other way than that specified in this manual.
- DO NOT use the device if you have had a recent surgery or other surgical procedures.
- DO NOT use the device over a bacterial infection as this may result in the worsening of symptoms.
- DO NOT use the device on painful areas. Consult your General Practitioner when you experience painful areas during use.
- DO NOT use the device on, near, or by children or those with reduced physical and mental abilities. Close supervision is necessary when this device is used, cleaned or stored near children or those with reduced physical or mental abilities.
- DO NOT use the device with any other USB-C cable the one provided by NOOANCE.
- DO NOT use the device with any other topical other than specified by NOOANCE.

Note: The device should only be supplied with 5V DC low voltage. If any part of the main cable is damaged, discontinue use. Contact customer support to arrange a cable replacement.

⚠ Precautions

- Do not scratch or damage the light-emitting surface.
- It will void your warranty if you do not operate the device in accordance with the instruction in the manual.
- Read all instructions carefully before using this device.
- Store the device according to the protocol worked out in this manual.
- This device is intended to be used by one person.
- The device is designed for cosmetic use only, for individuals with good health. Please contact your General Practitioner if you are unsure about your health status.
- The battery of the device is not removable by the user. If the device is to be disposed of, do not put into a fire. Dispose of the device per local regulations.
- Inspect the device and do not use the device if you see any signs of damage such as internal battery leaks. Contact our customer support.
- The device should only be supplied with 5V DC low voltage. If any part of the main cable is damaged, discontinue use. Contact customer support to arrange a cable replacement.
- Consult your General Practitioner before using the device if you have a fever, acute inflammation, a severe circulatory disorder, light sensitivity, anaesthetic, or desensitized areas.
- We highly suggest taking caution when taking any medications. If in doubt, consult your General Practitioner.

Travel And International Use

The device's lightweight design makes it a perfect travel companion. Here are some travel tips for you:

- If you are travelling to a different region of where you bought the device, please contact the retailer you purchased the LED Neck Mask from. The NOOANCE LED Neck Mask is only supplied with a 5V DC Low Voltage.
- Confirm with the aircrew that the device can be carried and/or used on the airplane prior to flight.

Disclaimer

DO NOT use the Device to treat any other conditions apart from the conditions that are indicated in the Indications for Use statement in Chapter General Device Information.

NOOANCE Led Mask Multi-Functions is not intended for use in the diagnosis of disease or other conditions or in the cure, mitigation, treatment or prevention of disease, nor is it intended to affect the structure or function of the body.

Risks

Potential Risks

- Immediately after treatment skin may appear red Skin should return to normal within 1-2 hours.
- Temporary headaches may appear after the light treatment. Consult you health professional if the headaches stay over a longer period of time.

Light Sensivity Test

We highly recommend users perform a light sensitivity test before their first use of the device.


1. Place the device on the inner forearm.
2. Turn on the device to start a treatment.
3. Wait till the treatment is done. The treatment will take 5 minutes to complete.
4. If your skin feels hot or uncomfortable during the treatment directly stop! You might be sensitive to light.
5. The device will automatically turn off after the treatment is done. Remove the device and wait for at least 24 hours.

Check your skin regularly over those 24 hours, if your skin appears red and/or itchy the device is unsuitable for you. In the case of light sensitivity please contact customer support.

How To Use The Device (Two ways to use)

Preparation: Before starting the treatment, please thoroughly cleanse your skin.

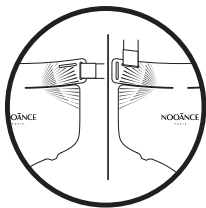Attach the adjustable straps to the device by threading through the sides of the device.

Step 4

If you wish, you can apply a moisturizing serum to your skin after the treatment.

Connect the controller to use

NECK

CHEST


HANDS

LEGS


LOWER ABDOMEN

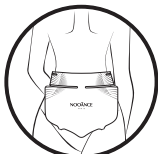BACK


Step 1

Apply the device to the desired body part by threading the straps through the appropriate sides.

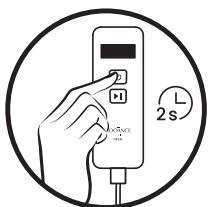Step 2

Long press the ON/OFF button for 2 seconds to turn on the device on the Controller.

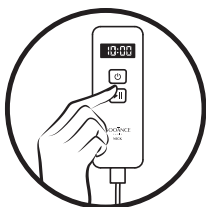Step 3

Press the start/pause button to begin the 10-minute treatment. A "beep" indicates the treatment has started. Once the treatment is complete, the device will turn off automatically with two "beep" sounds.

Connect the adapter to use

NECK

CHEST

HANDS

LEGS


LOWER ABDOMEN

BACK

Step 1

Apply the device to the desired body part by threading the straps through the appropriate sides.

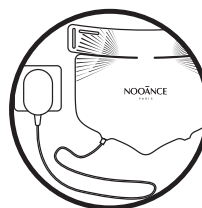Step 2

Connect the device to the adapter to begin the 10-minute treatment. (Manual timing is required.)


Step 3

If you wish, you can apply a moisturizing serum to your skin after the treatment.

Installation of Straps

There are 3 straps with 2 sizes, one is 25cm, the other two are 80cm. Attach the adjustable straps to the device by threading through the sides of the device.

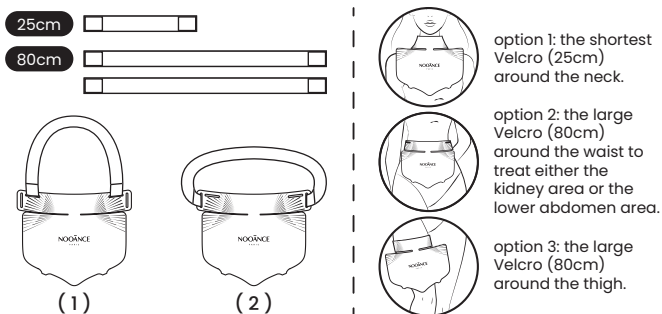Treatment Protocol

Use the device 1 time per day for a minimum of 4 weeks. After that you can continue using the device as required based upon your lifestyle and skin condition. It is not recommended to use the device more than 7 times per week.

Charging

Charge your NOOANCE device fully for 4 hours before your first use. Only use a 5V adapter or computer USB port (avoid fast chargers as they can damage the battery).

To charge the device insert the USB-C cable into the controller. Insert the USB-A side of the cable into a USB power supply. When the controller is charged, the LED indicator on the device flashes white, and when the controller is fully charged, the indicator will be steady white.

Before charging, make sure that the plug and socket are completely dry. Failure to do so may result in electric shock, short circuit, or fire.

Note: Always unplug the device before use, the device will not work while being charged.

Care And Storage

- The lifetime of the Lithium battery can be maximized by allowing it to fully charge and then fully discharge during regular use.
- Store the device at room temperature in dry conditions.
- Store the device in the supplied packaging. Make sure the device will not bend, stretch or crush while being stored or when travelled with.
- Avoid storing the device in dusty conditions, in direct sunlight or in damp and humid conditions.
- This device contains no serviceable parts.
- Before each use, with the device switched off, clean the device following the cleaning instructions.

Cleaning

- Clean the device with a water-based wipe or damp clean cloth. Stubborn stains can be removed by wiping gently with 70% alcohol.
- DO NOT use any other cleansing devices as these may damage the device.
- DO NOT clean the device while it is turned on.
- DO NOT clean the device while charging.

Warranty

NOOANCE Led Mask Multi-Functions is warranted to be free from defects in material and workmanship for 24 months from its original purchase date under normal use. This warranty extends only to the original retail purchaser with original proof of purchase and only when purchased from an authorized retailer or reseller. This warranty only extends if the device is used in conjunction with authorized device accessories. If the device should become defective within the warranty period, contact the retailer that you purchased the device from. This warranty does not cover devices damaged by:

- Accident, misuse, abuse, or alteration.
- Use with unauthorized accessories or use other than as instructed.
- Connecting to incorrect current and voltage.

Troubleshooting

If the actions suggested fail to fix the problem, or the problem reoccurs, and your device is within the warranty period please return to Nooance or contact to by email bonjour@nooance-paris.com.

- ⚠ DO NOT attempt to investigate any further.
- ⚠ DO NOT attempt to modify the device in any way. The device has NO parts you can repair yourself.

Problem	Cause	Solution
The device doesn't switch on.	The device is not fully charged, because the charging cable was not properly inserted into the power outlet.	Insert the charging cable properly into the power outlet.
	The device is not fully charged, because the charging cable was not properly inserted into the device.	Insert the charging cable properly into the device and charge the device.
	The power outlet didn't work during charging.	Try a different 5V adapter or power outlet.
	After pressing the On/Off for a full 2 seconds, there is no short press of the "Start/Pause" button for 1 second.	After pressing on/off for a full 2 seconds, you need to press the "Start/Pause" button for 1 second to start treatment.
The device switches off unexpectedly.	The battery level is too low.	Charge the device.

⚠ Disposal

DO NOT dispose of the device with household waste. The Device should be recycled. The device contains a Lithium-ion battery. For protection of the environment, please safely dispose of the device following your local guidelines.

Symbols

Symbol	Explanation
	Indicates the device is CE Marked and authorized to sell in Europe. (Lower-Voltage-Directive)
	The Waste of Electrical and Electronic Equipment (WEEE), Indicates the device may not be discarded as household waste (2012/19/EU).
	Indicates that the electromagnetic radiation from the device is below the limits specified by the Federal Communications Commission and the manufacturer has followed the requirements of the Supplier's Declaration of Conformity authorization procedures.
	This warning symbol on the device label is a warning to read the enclosed documents before use and take all safety precautions required.
	Consult instructions for use.
	Date of Manufacture.
	This symbol means keep the device dry.
	This symbol indicates the serial number of the device.
	This symbol refers to the batch code of the device.
	This symbol indicates the USB connector.
	This symbol means that the device uses a Lithium-ion Polymer battery.
	This symbol means keep the device away from direct sunlight.
	This symbol means the upper and lower operating temperatures that the device can be safely used at.
	This symbol indicates the upper and lower operating humidity that the device can be safely used at.

	This symbol means that the cardboard packaging that the device came in can be recycled.
	Type BF Applied Part, B indicates it is used for the body and F indicates that it is not connected to the mains voltage.
	Unique Device Identifier.
	This symbols implies conformity to UK Directives on health, safety and environmental protection standards. UKCA stands for UK Conformity Assessed.
	Indicates the device is CA65 Marked and authorized to sell in the United States.
	Indicates the device is ROHS Marked and authorized to sell in Europe.


 Made in China

Designed in France


 NOOANCE INNOVATIONS
22 Rue Beaujon, 75008 Paris, France

**COU
NECK**

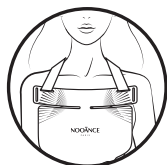
**TORSE
CHEST**

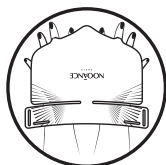
**MAINS
HANDS**

**JAMBES
LEGS**

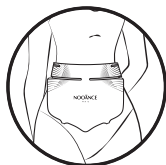**ABDOMEN
LOWER ABDOMEN**

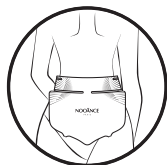**DOS
BACK**