NOOANCE

PARIS

MASQUE LED BEAUTÉ
BEAUTY LED MASK

ACCESS +

Français	—————	01-11
English	—————	12-22
Deutsch	—————	23-33
Español	—————	34-44
Italiano	—————	45-55
العربية	—————	56-66

Sommaire

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Table of Contents

The science of light therapy	13
Contents of the box	14
Usage	15
Before first use	
Frequency of use	
User Guide	
Product Information	17
Device Specifications	
Warnings	
Contraindications	
Warranty	
Air Travel	
Disclaimer	
Troubleshooting	
Disposal	
Symbols	21

The science of light therapy

The skin has the remarkable ability to absorb light energy. It uses this energy to boost the cells to produce more energy (ATP), which increases their performance and allows them to carry out their natural functions, but at a faster rate.

Cells supercharged by light therapy use this energy to boost the skin's fibroblasts to produce more collagen and elastin proteins and to stimulate the natural processes of cell renewal.

The benefits of light therapy :

- boost cell regeneration
- stimulate skin metabolism
- repair and soothe irritations

The NOOANCE LED Mask features 198 LEDs with a precision of ± 2 nanometers :

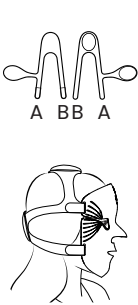
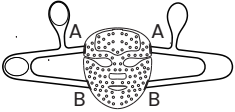
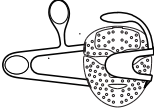
- 99 red LEDs with a wavelength of 633 nanometers
- 99 near-infrared LEDs with a wavelength of 830 nanometers.

Red and near-infrared LEDs help to :

- reduce wrinkles
- improve firmness
- give an immediate radiance boost

A powerful synergy of 633 nm and 830 nm wavelengths delivering visible results to revitalize cells and improve your skin's radiance, firmness and health.

Contents of the box

		Additional information
<p>NOOANCE LED ACCESS +</p>		<ul style="list-style-type: none"> ● Clean the device only with a clean, damp cloth or water wipe. You can also use an antibacterial cleaner or wipe. ● DO NOT clean the appliance when it is switched on or charging.
<p>Adjustable Straps</p>		<ul style="list-style-type: none"> ● Installation of straps <ul style="list-style-type: none"> ○ Sticky surface  <p style="text-align: center;">①</p>   <p style="text-align: center;">② ③</p> <ul style="list-style-type: none"> ● Straps marked (A) are attached to the top slots and straps marked (B) to the lower slots.
<p>USB cable</p>		<ul style="list-style-type: none"> ● Use the USB cable to charge the mask, the indicator light flashes during charging. Note that the mask does not function while charging.

Silicone Mask Protectors		<ul style="list-style-type: none"> ● When using a LED mask, red light and infrared light pose no danger to the eyes. Silicone mask protectors are optional but can be worn for added comfort. ● The red light of the LEDs may appear orange after a few seconds of use as the eye adapts to this brightness. This is completely normal.
Manual		<ul style="list-style-type: none"> ● Read this manual carefully before using your mask.

Usage

Before first use :

We strongly recommend that users carry out a light sensitivity test before first use.

1. Place the device on the inside of the forearm.
2. Switch on the device to start your session.
3. Turn the mask off after 5 minutes (the timer is set for 10 minutes by default).

If your skin becomes excessively hot or itchy during the session, stop immediately. You may be sensitive to light.

Check your skin regularly for the next 24 hours. If your skin is red and/or itchy, the appliance is not suitable for you. Please contact customer service.

Frequency of use

We recommend using the mask for 10 minutes per day, between 5 and 6 times per week.

User Guide

Preparation: Before starting the treatment, please thoroughly cleanse your skin.

Usage: Apply the mask to your face, and strap it around your head using the adjustable straps.

Tip: To optimise the results of your Nooance mask, we recommend that you use one of our aqueous serums recommended on our website www.nooance-paris.com

At the end of the session, you can also apply a moisturising cream or serum.

Please insert the protective glasses for more comfort.

Press 3 seconds the button at the top of the mask to start the 10-minute session.

A beep indicates that the session has begun.

Once your session is over, the mask will automatically turn off with two "beeps".

Product Information

Device Specifications

LED wavelengths	633 nm \pm 2 nm (Red) 830 \pm 2 nm (Infrared)
Max. output optical power density	45mW/cm ²
Energy delivered	27J/cm ²
Operating temperature	0°C to 40°C 32°F to 104°F
Storing and transport temperature	-5°C to 55°C 23°F to 131°F
Operating humidity	30% to 95% without condensation
Storage and Transport Humidity	20% to 95% without condensation
Charging supply	5V \equiv 1A
Charging port	USB-C
Charging cable	USB-C to USB-A
Charging cable length	1.2 meters / 47.244 feet
Charging adapter power output	1A
Battery	3.7V, 600mAh
Battery type	Lithium-ion
Charging time	90 minutes / 1.5 hours
Recommended use	50 to 60 minutes per week
Timer	10 minutes
Dimensions of the device	180.33*213.61*104.07 mm
Weight of the device	0.275 kg / 9.7 oz

Warnings

- Read all the instructions carefully before using this device.
- Stop using the device and consult your doctor if you experience any adverse reactions such as discomfort, dizziness, skin discolouration, skin irritation, or any other reaction that seems abnormal.
- Stop using the device if it overheats or if you feel unpleasant heat on your skin.
- DO NOT allow children or people with reduced ability to play with the appliance. The appliance must not be cleaned or charged by children.
- DO NOT point the LED light directly into your eyes or those of another person.
- DO NOT use the appliance if it is not working properly, if it has been dropped, or if it has been damaged. Contact customer service for warranty and repairs.
- DO NOT use, charge, or store the appliance near hot surfaces, in contact with water, or in direct sunlight.
- DO NOT attempt to repair the device. Repairing or opening the appliance will invalidate your warranty. The appliance contains no serviceable parts.
- DO NOT use the appliance with accessories not authorised by the manufacturer.
- DO NOT use solvents or abrasive cleaners on the appliance.
- DO NOT use the device in any context other than that specified in this manual.
- In very rare cases, temporary headaches may occur after a light therapy session. Consult your doctor if symptoms persist and discontinue use of the device.
- DO NOT scratch or damage the light-emitting surface.
- Inspect the device and do not use it if it shows signs of damage such as battery leaks. Contact our customer service.

- Your warranty will be void if you do not use the device according to the instructions in the manual.

Contraindications

- DO NOT use the device if you suffer from a photosensitive disorder (sensitivity to light), are taking photosensitising medication, or if you suffer from epilepsy. If in doubt, consult your doctor.
- DO NOT use the device if you have active cancer, a cancerous lesion, an open wound, warts, or any skin condition. Consult your doctor if in doubt.
- DO NOT use the device on a bacterial infection, as this may aggravate the symptoms.
- DO NOT use the device if you are pregnant, planning to become pregnant, or breastfeeding.
- DO NOT use the device on anyone under the age of 18.
- DO NOT use the device if you have a known allergic reaction to silicone. The device contains plastics and silicone that come into contact with the skin during use. If in doubt about allergies, contact your doctor before using the device.

Warranty

- The LED mask comes with a 24-month guarantee for the original purchaser, subject to proof of purchase from an authorised dealer. This guarantee covers normal use with approved accessories. In the event of a fault, contact our customer service department.
- Please note that the guarantee does not cover damage caused by accidents, incorrect use, use of non-recommended accessories, or poor electrical connection.

Air Travel

- The device is allowed in the cabin. If you wish to check it in, please consult your airline regarding Lithium batteries.

Disclaimer

- The NOOANCE Face Mask is not intended for use in the diagnosis of diseases and other conditions, nor for the cure, mitigation, treatment, or prevention of diseases, and is not intended to alter the structure or function of the body.

Troubleshooting

If the actions suggested fail to fix the problem, or the problem reoccurs, and your device is within the warranty period, please contact our customer service by writing to bonjour@nooance-paris.com.

⚠ DO NOT attempt to investigate any further.

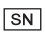
⚠ DO NOT attempt to modify the device in any way. The device has NO parts you can repair yourself.

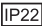
Problem	Cause	Solution
The device doesn't switch on.	The power outlet didn't work during charging.	Try a different 5V adapter or power outlet.
	The on/off button has not been pushed for a full 3 seconds.	Push the on/off button for 3 seconds straight.
The device switches off unexpectedly.	The battery level is too low.	Charge the Battery.

⚠ Disposal

DO NOT dispose of the device with household waste. The device should be recycled. The device contains a Lithium-ion battery. For protection of the environment, please safely dispose of the device following your local guidelines.

Symbols

Symbol	Explanation
	Indicates the device is CE Marked and authorized to sell in Europe. (Lower-Voltage-Directive)
	The Waste of Electrical and Electronic Equipment (WEEE), Indicates the device may not be discarded as household waste (2012/19/EU).
	Indicates that the electromagnetic radiation from the device is below the limits specified by the Federal Communications Commission and the manufacturer has followed the requirements of the Supplier's Declaration of Conformity authorization procedures.
	This warning symbol on the device label is a warning to read the enclosed documents before use and take all safety precautions required.
	Consult instructions for use.
	Date of Manufacture.
	This symbol means keep the device dry.
	This symbol indicates the serial number of the device.
	This symbol refers to the batch code of the device.
	This symbol indicates the USB connector.
	This symbol means that the device uses a Lithium-ion Polymer battery.
	This symbol means keep the device away from direct sunlight.
	This symbol means the upper and lower operating temperatures that the device can be safely used at.

	This symbol indicates the upper and lower operating humidity that the device can be safely used at.
	This symbol means that the cardboard packaging that the device came in can be recycled.
	IP22
	Type BF Applied Part, B indicates it is used for the body and F indicates that it is not connected to the mains voltage.
	Unique Device Identifier.
	This symbol implies conformity to UK Directives on health, safety and environmental protection standards. UKCA stands for UK Conformity Assessed.
	Indicates the device is CA65 Marked and authorized to sell in the United States.
	Indicates the device is ROHS Marked and authorized to sell in Europe.

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